

Document made available under the Patent Cooperation Treaty (PCT)

International application number: PCT/US04/041801

International filing date: 10 December 2004 (10.12.2004)

Document type: Certified copy of priority document

Document details: Country/Office: US
Number: 60/531,379
Filing date: 19 December 2003 (19.12.2003)

Date of receipt at the International Bureau: 13 January 2005 (13.01.2005)

Remark: Priority document submitted or transmitted to the International Bureau in compliance with Rule 17.1(a) or (b)



World Intellectual Property Organization (WIPO) - Geneva, Switzerland
Organisation Mondiale de la Propriété Intellectuelle (OMPI) - Genève, Suisse



THE UNITED STATES OF AMERICA

TO ALL TO WHOM THESE PRESENTS SHALL COME:

UNITED STATES DEPARTMENT OF COMMERCE

United States Patent and Trademark Office

January 05, 2005

THIS IS TO CERTIFY THAT ANNEXED HERETO IS A TRUE COPY FROM THE RECORDS OF THE UNITED STATES PATENT AND TRADEMARK OFFICE OF THOSE PAPERS OF THE BELOW IDENTIFIED PATENT APPLICATION THAT MET THE REQUIREMENTS TO BE GRANTED A FILING DATE.

APPLICATION NUMBER: 60/531,379

FILING DATE: *December 19, 2003*

RELATED PCT APPLICATION NUMBER: *PCT/US04/41801*



Certified By

Jon W Dudas

Under Secretary
of Commerce for Intellectual Property
and Acting Director of the
United States Patent and Trademark Office



121903

16563 U.S. PTO

Mail Stop Provisional Patent Application
 Commissioner for Patents
 PO Box 1450
 Alexandria, VA 22313-1450

Docket No: AM101242 L1
 Patent

PROVISIONAL APPLICATION COVER SHEET

Transmitted herewith for filing under 37 CFR 1.53(c) is the provisional patent application of the following Inventor(s): Ene ALLIKMETS
 For: SERUM-FREE VERO CELL BANKING PROCESS

22151 U.S. PTO
 60/531379



121903

1. Papers enclosed which are required for filing date under 37 CFR 1.51(c) and 1.53 (c):

- ☒ Pages of specification – 13 pages
☐ Sequence Listing – pages
☒ Pages of claims – 3 pages
☒ Page(s) of abstract – 1 pages
 Sheet(s) of drawing – 5 sheets
☐ Formal
☒ Informal

2. Inventor(s) name(s) and address(es) (or list attorney docket number as alphanumeric identifier):

1. Ene Allikmets 200 Maple Road, Cornwall on Hudson, NY 12520
 2.
 3.
 4.
 5.

3. Assignment

An assignment of the invention to:
 Wyeth
 Five Giralda Farms
 Madison, New Jersey 07940

- ☐ is attached under separate Recordation Form Cover Sheet.
☒ will follow.

CERTIFICATE OF MAILING 37 CFR §1.10

I hereby certify that this paper and the documents referred to as enclosed therein are being deposited with the United States Postal Service on the date written below in an envelope as "Express Mail Post Office to Addressee" Mailing Label Number EU238168010US addressed to the Mail Stop Provisional Patent Application, Commissioner for Patents, PO Box 1450, Alexandria, VA 22313-1450.

Date


12/19/03

RoseAnn Jastrow

4. Filing Fee \$160.00
5. Method of Payment of Fees:
Charge Deposit Account No. 01-1425 in the amount of \$160.00.
A duplicate of this transmittal is attached.
6. Instructions as to Overpayment/Underpayment:
Credit any overpayment and charge any underpayment to Deposit Account No. 01-1425.
7. Authorization to Charge Additional Fees
☒ The Commissioner is hereby authorized to charge the following additional fees by this paper and during the entire pendency of this application to Deposit Account No. 01-1425:
☒ 37 CFR 1.51(c)(4) and 1.53 (c) filing fees
☒ Surcharge for filing the basic filing fee and/or complete cover sheet on a date later than the filing date of the application.
8. SEND CORRESPONDENCE TO:

Customer Number: 25291
- DIRECT ALL TELEPHONE CALLS TO:

Name: Carol E. Rozek
Tel. No. (845) 602-4760
9. The invention was made by an agency of the U.S. Government or under a contract with an agency of the U.S. Government.
☒ No
☐ Yes, the name of the U.S. Government agency and the Government contract numbers are _____.


Carol E. Rozek
Attorney for Applicants
Reg. No. 36,993

Wyeth
Patent Law Department
Five Giralda Farms
Madison, NJ 07940
Tel. No. (845) 602-4760

SERUM-FREE VERO CELL BANKING PROCESS

BACKGROUND OF THE INVENTION

5 The present invention relates to a serum-free cell-freezing medium and the use of this medium in a process for generating stable serum-free cell banks for viral vaccine production.

Fetal bovine serum (FBS) is a cryostabilizing agent commonly used in cell banking. Serum-containing culture systems, however, are becoming undesirable for the large-scale production of vaccines. There are a number of disadvantages of serum supplementation, including batch-to-batch variation in composition, the high protein content that hinders product purification, and the potential for viral, mycoplasma, or prion contamination. Furthermore, the recent threat to human health caused by the undefined agents of bovine spongiform encephalopathy (BSE) is likely to limit the continued use of bovine serum in culture processes used for the synthesis of health care products such as viral vaccines [Butler, et al., *Biotechnol. Prog.*, 16, 854-858 (2000)]. Therefore, a cryostabilizing agent of non-animal origin is needed to replace FBS.

20 SUMMARY OF THE INVENTION

The present invention provides a serum-free cell-freezing medium consisting essentially of a virus-production serum-free medium (VP-SFM) supplemented with (a) an enzymatic hydrolysate cryostabilizer selected from the group consisting of soy hydrolysate and rice hydrolysate, added at 2 g per liter of said medium, and (b) dimethylsulfoxide (DMSO). In one embodiment, the DMSO is added at 10 g per liter of said medium.

The present invention further provides for the use of this serum-free cell-freezing medium in a process that generates standardized and consistent master and working cell banks, which are used as raw material in vaccine production. This process comprises the steps of:

(a) initiating a culture, which comprises thawing frozen Vero cells and adding them to growth medium in a T-150 cm² flask, wherein the growth medium consists of VP-SFM with 4mM L-glutamine, incubating the cells overnight at 37°C and 5% CO₂, and refeeding the culture with fresh growth medium the following day;

- (b) propagating and amplifying the cells, which comprises growing the cells to confluence in the T-150 cm² flask incubated at 37°C and 5% CO₂, removing the medium, washing the flask phosphate buffered saline (PBS) without calcium and magnesium, adding trypsin to the flask and incubating at room temperature for a
- 5 sufficient time to dislodge the cells from the flask, neutralizing the trypsin with soybean trypsin inhibitor (STI) and adding VP-SFM for nutritional support, seeding the resulting suspension into five T-150 cm² flasks and adding VP-SFM to each flask to a level of 50 ml, incubating the cells at 37°C and 5% CO₂ for three to four days, pooling the cell suspensions from the five flasks, seeding a cell factory with the
- 10 pooled suspension, refeeding the cell factory, and harvesting the cell factory; and
- (c) freezing the cell bank, which comprises centrifuging the harvested cells from the cell factory for 10 minutes at 210xg at 4°C, resuspending the cells in the serum-free cell-freezing medium of claim 1 at a density of 2×10^6 to 2×10^7 cells/ml, dispensing the cell suspension into cryovials at one ml of cell suspension per vial,
- 15 freezing the cells using an active rate control freezer, and storing the cells in liquid nitrogen,
- wherein the stable serum-free Vero cell bank thus produced has a cell viability of at least 80% and a recovery doubling time between 40 and 60 hours after one year.

20 In a specific embodiment, this process comprises the steps of:

- (a) initiating a culture, which comprises thawing 2×10^7 frozen Vero cells and adding them to 50 ml of growth medium in a T-150 cm² flask to obtain a cell density of $4\text{-}5 \times 10^5$ cells/ml, wherein the growth medium consists of VP-SFM with 4mM L-glutamine, incubating the cells overnight at 37°C under 5% CO₂, and refeeding the
- 25 culture with fresh growth medium the following day;
- (b) propagating and amplifying the cells, which comprises growing the cells to confluence in the T-150 cm² flask incubated at 37°C under 5% CO₂, removing the medium, washing the flask two times with 20 ml phosphate buffered saline (PBS) without calcium and magnesium, adding 5 ml trypsin to the flask and incubating at
- 30 room temperature for a sufficient time to dislodge the cells from the flask, neutralizing the trypsin with 5 ml of soybean trypsin inhibitor (STI) and adding 10 ml of VP-SFM for nutritional support, seeding the resulting suspension into five T-150 cm² flasks at a concentration of 4×10^4 cells/cm² and adding VP-SFM to each flask to a level of 50

ml, incubating the cells at 37°C under 5% CO₂ for three to four days, pooling the cell suspensions from the five flasks, seeding a cell factory with the pooled suspension, refeeding the cell factory, and harvesting the cell factory; and

- (c) freezing the cell bank, which comprises centrifuging the harvested cells from the cell factory for 10 minutes at 210xg at 4°C, resuspending the cells in the serum-free cell-freezing medium described above at a density of $1-2 \times 10^7$ cells/ml, dispensing the cell suspension into cryovials vial at one ml of cell suspension per vial, freezing the cells using an active rate control freezer, and then storing the cells in liquid nitrogen,
- wherein the stable serum-free Vero cell bank thus produced has a cell viability of at least 80% and a recovery doubling time between 40 and 60 hours after one year.

The present invention further provides a stable serum-free Vero cell bank for vaccine production produced by the process immediately described above, wherein the average cell bank viability is at least 80% and the average recovery doubling time is between 40 and 60 hours after one year.

BRIEF DESCRIPTION OF THE DRAWINGS

FIG. 1 is a scale-up scheme for 100 vial cell banks. Cell banks of 200 vials and 400 vials follow the same scheme, using two cell factories and four cell factories, respectively.

FIG. 2 is a graph comparing the cell viabilities of 10 vial serum-free cell banks containing plant-derived cryostabilizers with that of FBS-containing cell banks.

FIG. 3 is a graph comparing the recovery doubling times of 10 vial serum-free cell banks containing plant-derived cryostabilizers with that of FBS-containing cell banks.

FIG. 4 is a graph comparing the cell viabilities of large-scale cell banks (100, 200 and 400 vials), which were made using Hy-Soy® and HyPep® Rice as the cryostabilizers.

FIG. 5 is a graph comparing the recovery doubling times of large-scale cell banks (100, 200 and 400 vials), which were made using Hy-Soy and HyPep Rice as the cryostabilizers.

DETAILED DESCRIPTION OF THE INVENTION

Implementation of a correct and standardized cell banking process is fundamental to the provision of cell stocks with consistent performance and characteristics. Cell banking should include a master cell bank of homogeneous source cells and numerous working cell banks generated from the master cell bank. The number of vials of a given bank will vary depending on the intended use, but a typical master cell bank may contain 200 vials and a typical working cell bank may contain up to 400 vials, utilizing the banking process of the present invention. This process ensures a long-term supply of standardized cells that will provide a stable substrate for consistent production.

Vero cells (immortalized African Green Monkey kidney cells, ATCC # CCI-81) have been known and used in the art for the production of human viral vaccines for many years. The cells are well characterized and have an excellent safety record.

The serum-free Vero cell banking process of the present invention consists of three major steps: initiating the culture, propagating and amplifying the cells, and freezing the bank. Optimization of each of these steps has made it possible to generate high quality serum-free cell banks. "High quality" means the viability of the cells is high and the recovery is fast, i.e., the recovery doubling time is low. The combination of these cell growth and freezing conditions results in improved recovery of Vero cells. Two of the most important recovery parameters for consistent production of cells from a cell bank are the recovery time and the thaw viability. The recovery time is closely related to the plating efficiency but it also provides more information about the performance of the cells after thawing. Both of these parameters were measured and used to optimize the instant cell banking process.

Phenol red was eliminated from the cell-freezing medium due to reports that this component could have estrogenic effects on the cells. [Berthois, et al. *Cell Biology*, 83, 2496-5000 (1986)]. The composition of the cell-freezing medium was optimized to contain 2 g/L Hy-Soy (or HyPep Rice), 10 g/L DMSO, and the remainder VP-SFM as a growth medium.

As used herein "VP-SFM" means Gibco™ VP-SFM (commercially available from Invitrogen Corp., Carlsbad, CA) modified to contain the ingredients listed in Table 1.

Table 1**GIBCO VP-SFM
(modified)**

5

L-Alanine	Sodium Phosphate Dibasic Anhyd	Putrescine 2HCl
Glycine	Magnesium Sulfate Anhyd	Recombinant Human Insulin
L-Arginine HCl	Calcium Chloride Anhyd	Epidermal Growth Factor (EGF), Recombinant
L-Asparagine H ₂ O	D Calcium Pantothenate	Sodium Phosphate Dibasic Anhyd.
L-Aspartic Acid	Magnesium Chloride Anhyd	Sodium Phosphate Monobasic H ₂ O
L-Cysteine HCl H ₂ O	Ascorbic Acid 2 Phosphate Magnesium Salt	HEPES
L-Cystine 2HCl	Sodium Chloride	Ferric Citrate Complex
L-Glutamic Acid	Vitamin A Acetate	Adenine Sulfate
Glutathione Reduced	Vitamin D ₂	Adenosine-5-Triphosphate
L-Histidine HCl H ₂ O	Menadione	D-Glucose (Dextrose)
L-Isoleucine	DL Lipoic Acid Thioctic	Vitamin B ₁₂
L-Leucine	Linoleic Acid	Choline Chloride
L-Lysine HCl	Uracil	I-Inositol
L-Methionine	Thymidine	Niacinamide
L-Phenylalanine	Biotin	Pyridoxal HCl
L-Proline	Folic Acid	2-Deoxyribose
L-Serine	Riboflavin	Ethanolamine HCl
L-Threonine	Sodium Hypoxanthine	Phosphoethanolamine
L-Tryptophan	Sodium Xanthine	Potassium Chloride
L-Valine	Para Amino Benzoic Acid	Sodium Pyruvate
Thiamine HCl	Nicotinic Acid	Proprietary plant hydrolysate
L-Tyrosine Disodium Salt	Pyridoxine HCl	Yeast Extract UF

As a process improvement, the culture was initiated from double the cell density that has been customary in the art. Therefore, rather than starting with one vial, the culture is initiated from two frozen vials from the source bank thawed into one T-150 cm² flask. Warm VP-SFM is added drop wise to avoid shock or damage to the cells. Refeeding the culture the next day after thawing removes remaining cryostabilizer and cell debris, and improves the recovery.

The scale-up scheme and schedule of the present invention enables the propagation of Vero cells with a minimal number of passages (e.g., no more than four passages) and manipulations, while still keeping the cell concentration consistently high to maintain a healthy culture. The scale-up scheme is: two frozen vials from the

source bank thawed into one T-150 cm² flask, expanded into five T-150 cm² flasks that are used to inoculate one 10-tray cell factory. A 100-vial cell bank can be generated from one cell factory. This scheme allows for easy adjustments in bank size.

5 The trypsinization scheme was also optimized. Two buffer rinses of the cell monolayer were found to be optimal for good recovery after thawing. Gamma-irradiated porcine trypsin was chosen to minimize possible viral contamination. Because sufficient heat transfer to all ten levels of the cell factory is questionable, creating unnecessary microenvironments, the trypsinization is performed at room
10 temperature with chilled trypsin.

 The combination of freezing process conditions improved recovery of the bank considerably. First, the harvested and pooled cell suspension is centrifuged at 210xg for 10 minutes as opposed to 5 or 7 minutes, which causes considerable cell loss. Resuspension is performed gently with cold freezing medium. Second,
15 dispensing the cells into cryovials with a repeater pipettor made the dispensing procedure faster and reduced the time cells had to be exposed to the freezing medium before the actual freezing cycle. And third, the cell suspension during dispensing into cryovials and the already filled vials were cooled down to 2-8°C.

 Two aspects of the cryopreservation of Vero cells needed to be addressed for
20 the banking process: the optimum cell density in the cryovial, and programming the rate control freezer for good recovery of the cells upon thawing. During the phase change, as a material freezes the heat of transition from liquid to solid is released (the heat of fusion). The program was optimized so that the heat of fusion would be accounted for, forcing the sample to freeze at a constant rate. It was clearly shown
25 that the optimum cell density was $1-2 \times 10^7$ cells/ml and that the modifications to the freezing program were contributing to the good recovery of cells.

 The serum-free Vero cell banking process of the present invention provides a standardized and consistent way to generate reliable and stable cell banks for viral vaccine production. Animal-derived products were replaced in the growth and
30 freezing media by plant-derived substances such as Hy-Soy and HyPep Rice, which serve as cryostabilizers. The results show that serum-free cell banks using either 2% Hy-Soy or 2% HyPep Rice are stable for at least one year.

The improvements resulting in the instant cell banking process increase the viability of the cells upon thawing and reduce recovery time, which allows a more precise schedule for manufacturing and more consistent processes. It is also now possible for the manufacturing process to start from a smaller number of frozen vials.

- 5 Therefore one full-size 400-vial bank will last longer and consecutive runs will be more consistent.

Having now generally described the invention, the same will be more readily understood through reference to the following examples that are provided by way of illustration, and are not intended to be limiting of the present invention, unless
10 specified.

EXAMPLES

Several agents were tested as cryostabilizers in this study. Cell banks
15 containing Hy-Soy, HyPep Rice, HyPep Wheat and methylcellulose were compared with FBS-containing cell banks, and cell bank stability was examined. It was demonstrated that Hy-Soy and HyPep Rice were as effective as FBS as cryostabilizers for the first year of the study. Therefore, Hy-Soy can be used in serum-free cell banks and HyPep Rice can be used as an alternative.
20

Materials and Methods

Personal Protective Equipment: Sterile gloves, sterile sleeves, sterile jumpsuit, clean
25 room head covers, clean room beard covers, clean room foot covers.

Disposable Equipment: Disposable serological pipettes; 1.8 ml Sarstedt cryovials; Corning Costar 75 cm² and 150 cm² tissue culture flasks; ten tray Nunc cell factories.

30 *Glassware and Accessories:* Glass bell with C-flex tubing for dispensing medium from 10 L and 20 L medium totes; Gelman filters with C-flex tubing; cell factory adapters; 500 ml, 1 L, 2 L and 5 L bottles with transfer apparatus and tubing; hand bulb.

Reagents: Gibco™ VP-SFM formulated to contain the ingredients in Table 1; Gibco soybean trypsin inhibitor (STI); Gibco L-glutamine; Gibco phosphate buffered saline (PBS) without calcium and magnesium; Gibco bovine serum albumin (BSA); Hy-Soy, HyPep Rice, and HyPep Wheat from Quest International (Norwich, NY); Sigma (St. Louis, MO) dimethylsulfoxide (DMSO) and methylcellulose; 70% isopropyl alcohol; Gibco 0.25% irradiated trypsin/citrate; Gibco 0.025 M sodium citrate buffer; Hyclone (Logan, UT) Fetal Bovine Serum (FBS) or irradiated FBS (irFBS); Gibco Dulbecco's Minimal Essential Medium (DMEM); water for injection or deionized water; Gibco 0.05% Trypan Blue Solution.

Equipment: Bio safety hood; laminar flow hood or bio safety cabinet; spinner flask base; automatic (repeater) pipettor; microscope; inverted microscope; hemacytometer; circulating water bath; Forma Cryomed rate control freezer (RCF); peristaltic pump; dosing pump; Forma incubator at 37°C and 5% CO₂; Forma centrifuge; 216 swinging bucket type rotor.

Example 1

Thawing and Cell Growth

Two frozen stock vials containing 2×10^7 Vero cells were placed in a 37°C water bath. The vials were held in the water bath with agitation, and checked frequently until frozen cells were thawed. The outside of the vials were sanitized with 70% isopropyl alcohol and placed in a laminar flow hood. Thawed cells were quickly transferred to a T-150 cm² flask and 50 ml of growth medium was added drop wise while rocking the culture container. Growth medium consisted of VP-SFM with 4 mM L-glutamine. After approximately half of the medium was added, the rest was added at a slightly faster rate. The cells were placed in a Forma incubator (37°C, 5%CO₂) and allowed to attach overnight. The medium was aspirated off and 50 ml of fresh growth medium was added to the flask the following day.

Example 2

Cell Passing

Cells were grown to confluence (4 or 5 days) in the T-150 cm² flask incubated at 37°C and 5% CO₂ in the Forma incubator. The medium was aspirated off and the

flask was washed two times with 20 ml PBS, without magnesium and calcium. Five ml of trypsin was added to the flask and the flask was incubated at room temperature for 2 to 3 minutes to remove the cells from the flask. The trypsin was neutralized with 5 ml of STI and 10 ml of VP-SFM was added for nutritional support until the procedure was completed. The cells, now in a suspension, were sampled for counting. One ml of Trypan-blue solution was added to a one ml cell sample. This was mixed gently with a vortexor. Ten μ l of mixture was placed in the chamber of a hemacytometer. Cells resisting the stain were counted as live and the total cell count was determined.

The cell suspension was then seeded into five new T-150 cm² flasks at a concentration of 4×10^4 cells/cm². VP-SFM was added to each T-150 cm² flask to a level of 50 ml, and then the flasks were placed in the Forma incubator at 37°C and 5% CO₂. Three days later the flasks were refed. The flasks were incubated for a total of five days. The cell suspensions from the five flasks were pooled, and the pooled suspension was used to seed a Nunc cell factory. One confluent cell factory produces enough cells for a 100-vial cell bank.

Example 3

Cell Factory Set-up and Seeding

One of the adapters provided by Nunc for the cell factory was connected to a 10-inch piece of silicone tubing (ID 3/16", OD 3/8"). A tubing connector was used to attach an 8-inch length of C-flex tubing (ID .125", OD .200") to the silicone tubing connected to the cell factory. This is the inlet/outlet adapter. Both ends of tubing were wrapped with bioshield. The second adapter provided by Nunc for the cell factory was attached to a three-inch piece of silicone tubing (ID 3/16", OD 3/8") with a small Gelman filter (0.22 μ m). The adapter end was wrapped with bioshield. Both adapters were autoclaved for 30 minutes with a 15-minute dry time. In the BioSafety Cabinet (BSC), the adapters were attached to the cell factory as per the instruction manual provided by Nunc. Either port on the cell factory may be used for the inlet/outlet line and filter assembly line. A 10-L bag of VP-SFM or 3 x 1L bottles were obtained from the chill room and warmed to 37°C. Five L C-flex transfer lines were placed in sterile 5-L bottle for cell seeding. All bags used were equipped with C-flex by the manufacturer. The Sterile Connection Device (SCD) was used to connect the

medium to the seeding bottle and 2.5 L of medium was dispensed into the seeding bottle. Cell suspension was added to the seeding bottle. The SCD was used to connect the seeding bottle to the cell factory. The cell factory was placed in the fill and drain position. The medium with cells was slowly hand pumped into the cell
5 factory. The silicon tubing connected to the inlet/outlet line was immediately clamped off. A 0.22 μ m Gelman filter (Acro 50) with C-flex was connected to the cell factory using the SCD. The cell factory was placed in the level check position. When liquid was evenly distributed into the ten levels, the cell factory was placed into the horizontal position, and placed in a 37°C, 5% CO₂ incubator.

10

Example 4

Cell Factory Refeed

The cell factory was refed on days 2 and 5 after seeding. Using the SCD, the void bag was connected to the cell factory. The cell factory was gently placed in the
15 fill and drain position and liquid level marked. The void bag was placed on the floor to help drain the medium. It was necessary at times to pressurize the cell factory. A hand bulb was connected to the filter assembly to gently pressurize the cell factory. Once flow began, the hand bulb was disconnected. When medium was emptied, the bag containing the refeed medium was connected to the cell factory using the SCD
20 and the same volume of medium was slowly hand pumped into the cell factory. The cell factory was placed in the level check position. After level check was completed the cell factory was placed in its horizontal working position in the 37°C, 5% CO₂ incubator.

25

Example 5

Cell Factory Harvest

The cell factory was drained as described above. The first wash consisted of either 300-ml of 25mM sodium citrate buffer or 300 ml PBS without calcium or magnesium. The first wash was pumped into the cell factory. The wash was allowed
30 to equalize in each level. The cell factory was shaken gently to rinse all the cell surfaces. The first wash was drained and a second wash was applied in the same manner as the first wash. The second wash was drained and a bottle of trypsin solution (1.4 ml of Gibco 2.5% trypsin in 300 ml of 2.5 mM sodium citrate buffer) was

attached to the cell factory. The trypsin solution was pumped into the cell factory and the cell factory was gently rocked so that all the cell surfaces were covered. The cell surface was examined every 2 or 3 minutes to see if the cells had begun to detach. When cells began to detach from the cell factory surface, the cell factory was

5 rigorously tapped to release the rest of the cells. This was done until the cell factory appeared clear. The cell factory was then connected to a cell collection bottle containing 100 ml of 5 mg/ml soybean trypsin inhibitor (STI). The cell mixture was drained into this bottle. The medium wash bottle was connected to the cell factory, care being taken to clamp the cell collection bottle. The bottle containing 600 ml of

10 growth medium was attached to the cell factory to rinse the surfaces. The cell factory was reconnected to the cell collection bottle and the rinse medium was drained. This wash can contain up to 1/3 of the total cells.

Example 6

Freezing Protocol

15 Cells harvested from flasks or Nunc cell factories were sampled for counting, pooled, and centrifuged for 10 minutes at 210xg at 4°C in the Forma centrifuge. Cell-freezing medium was prepared consisting essentially of VP-SFM, supplemented with 10% DMSO, and 2% of one of the cryostabilizers tested in this study. The cells were

20 kept on ice during the process. The cells were then resuspended in the freezing medium at a density of 1.5×10^7 cells/ml. The cell concentrate was put into 1.8 ml Sarstedt cryovials at one ml of cell suspension per vial. Cells were frozen using an active rate control freezer. After freezing, the cells were stored in liquid nitrogen.

Example 7

Cell Banks

25 *Cell banks:* Initial experiments examining various cryostabilizers were accomplished using 10 vial cell banks derived from T-flasks. Based on the results of the 10 vial cell banks, the promising cryostabilizers were used in larger cell banks of

30 100, 200 and 400 vials. The cells needed for large-scale cell banks were grown in cell factories.

Cell expansion for large-scale cell banks: The expansion of cells for cell banks was performed as follows. Two vials of frozen cells were thawed for each 100

vials of the cell bank to be made. Two vials were used to inoculate one T-150 cm² flask. The T-150 cm² flask was refed the following day. After 4 days growth, the cells were split into five T-150 cm² flasks. Three days after the split, the flasks were refed. The flasks were incubated for five days. The flasks were harvested and used to
5 inoculate one Nunc 10 layer cell factory. The cell factory was refed after two and five days. On the seventh day, the cell factory was harvested. See FIG. 1 for the schedule of steps used to prepare the cell banks.

Cell bank testing: Five vials from each cell bank were thawed as described in Example 1. Each vial was used to seed two T-75 cm² flasks containing 25 ml of VP-
10 SFM. Cells were harvested from flasks when they were close to confluence and the following day. Doubling times were determined using cell count and time of growth. Cell viability was determined with Trypan-blue as described above.

RESULTS AND DISCUSSION

15 *Cell Bank Recovery:* Two parameters were used to determine cell bank stability: cell recovery doubling time (flasks harvested and counted when they were close to confluence and the following day) and the viability of cells after thawing. Cells were thawed immediately after freezing and at various time points. The cells were immediately tested for viability and placed in T-flasks. The T-flasks were
20 harvested to determine cell recovery doubling times.

Ten Vial Vero Cell Bank (See FIGS. 2 and 3): All cell banks had greater than 90% viability after being stored in liquid nitrogen for one year. The 10% DMSO only cell bank and the 0.1% methylcellulose cell bank dropped to below 80% at 15 months. It was also seen that the recovery of the methylcellulose bank was very
25 poor with doubling times of greater than 100 hours as compared to 45 to 50 hour doubling times of the other cell banks. The Hy-Soy, HyPep Rice and HyPep Wheat all had the doubling times and viabilities that were seen in the control cell banks (FBS- containing banks). Based on these data, large-scale cell banks of 100, 200 and 400 vials were made using Hy-Soy and HyPep Rice as cryostabilizers.

30 *Large-Scale Cell Banks* (See FIGS. 4 and 5): The one hundred vial Hy-Soy and HyPep Rice-containing cell banks' viability remained greater than 80% for nine months. Cell recovery also remained consistently between 45 to 60 hours doubling time. These results are consistent with the serum-containing cell banks (see FIGS. 2

and 3). The 200 vial Hy-Soy bank and the 400 vial Hy-Soy bank also showed a viability of over 80% and recovery doubling times between 40 hours and 60 hours over a twelve month period. The Hy-Soy and HyPep Rice large-scale cell banks are all comparable to the FBS-containing cell banks. From the above experimental
5 results, it can be concluded that serum-free cell banks using either 2% Hy-Soy or 2% HyPep Rice are stable for up to one year.

This invention has been described by a direct description and by examples. As noted above, the examples are meant to be only examples and not to limit the invention in any meaningful way. Additionally, one having ordinary skill in the art to
10 which this invention pertains in reviewing the specification and claims that follow, would appreciate that there are equivalents to those claimed aspects of the invention. The inventors intend to encompass those equivalents within the reasonable scope of the claimed invention.

What is claimed is:

1. A serum-free cell-freezing medium consisting essentially of a virus-production serum-free medium (VP-SFM) supplemented with (a) an enzymatic hydrolysate
5 cryostabilizer selected from the group consisting of soy hydrolysate and rice hydrolysate, added at 2 g per liter of said medium, and (b) dimethylsulfoxide (DMSO).
2. The serum-free cell-freezing medium of claim 1, wherein the DMSO is added
10 at 10 g per liter of said medium.
3. A process for generating a stable serum-free Vero cell bank, which process comprises the steps of:
 - (a) initiating a culture, which comprises thawing frozen Vero cells and adding
15 them to growth medium in a T-150 cm² flask, wherein the growth medium consists of VP-SFM with 4mM L-glutamine, incubating the cells overnight at 37°C and 5% CO₂, and refeeding the culture with fresh growth medium the following day;
 - (b) propagating and amplifying the cells, which comprises growing the cells to confluence in the T-150 cm² flask incubated at 37°C and 5% CO₂, removing the
20 medium, washing the flask phosphate buffered saline (PBS) without calcium and magnesium, adding trypsin to the flask and incubating at room temperature for a sufficient time to dislodge the cells from the flask, neutralizing the trypsin with soybean trypsin inhibitor (STI) and adding VP-SFM for nutritional support, seeding the resulting suspension into five T-150 cm² flasks and adding VP-SFM to each flask
25 to a level of 50 ml, incubating the cells at 37°C and 5% CO₂ for three to four days, pooling the cell suspensions from the five flasks, seeding a cell factory with the pooled suspension, refeeding the cell factory, and harvesting the cell factory; and
 - (c) freezing the cell bank, which comprises centrifuging the harvested cells from the cell factory for 10 minutes at 210xg at 4°C, resuspending the cells in the serum-
30 free cell-freezing medium of claim 1 at a density of 2×10^6 to 2×10^7 cells/ml, dispensing the cell suspension into cryovials at one ml of cell suspension per vial, freezing the cells using an active rate control freezer, and storing the cells in liquid nitrogen,

wherein the stable serum-free Vero cell bank thus produced has a cell viability of at least 80% and a recovery doubling time between 40 and 60 hours after one year.

4. A process for generating a stable serum-free Vero cell bank, which process
5 comprises the steps of:
- (a) initiating a culture, which comprises thawing 2×10^7 frozen Vero cells and adding them to 50 ml of growth medium in a T-150 cm² flask to obtain a cell density of $4\text{--}5 \times 10^5$ cells/ml, wherein the growth medium consists of VP-SFM with 4mM L-glutamine, incubating the cells overnight at 37°C and 5% CO₂, and refeeding the
10 culture with fresh growth medium the following day;
- (b) propagating and amplifying the cells, which comprises growing the cells to confluence in the T-150 cm² flask incubated at 37°C and 5% CO₂, removing the medium, washing the flask two times with 20 ml phosphate buffered saline (PBS) without calcium and magnesium, adding 5 ml trypsin to the flask and incubating at
15 room temperature for a sufficient time to dislodge the cells from the flask, neutralizing the trypsin with 5 ml of soybean trypsin inhibitor (STI) and adding 10 ml modified VP-SFM for nutritional support, seeding the resulting suspension into five T-150 cm² flasks at a concentration of 4×10^4 cells/cm² and adding VP-SFM to each flask to a level of 50 ml, incubating the cells at 37°C and 5% CO₂ for three to four days, pooling
20 the cell suspensions from the five flasks, seeding a cell factory with the pooled suspension, refeeding the cell factory, and harvesting the cell factory; and
- (c) freezing the cell bank, which comprises centrifuging the harvested cells from the cell factory at for 10 minutes at 210xg at 4°C, resuspending the cells in the serum-free cell-freezing medium of claim 2 at a density of $1\text{--}2 \times 10^7$ cells/ml,
25 dispensing the cell suspension into cryovials at one ml of cell suspension per vial, freezing the cells using an active rate control freezer, and storing the cells in liquid nitrogen,

wherein the stable serum-free Vero cell bank thus produced has a cell viability of at least 80% and a recovery doubling time between 40 and 60 hours after one year.

30

5. The serum-free cell-freezing medium of claim 1, wherein the enzymatic hydrolysate cryostabilizer is a soy hydrolysate.

6. The serum-free cell-freezing medium of claim 1, wherein the enzymatic hydrolysate cryostabilizer is a rice hydrolysate.

5 7. A stable serum-free Vero cell bank having a cell viability of least 80% and a recovery doubling time between 40 and 60 hours after one year.

8. A stable serum-free Vero cell bank having a cell viability of least 80% and a recovery doubling time between 40 and 60 hours after one year, wherein the cell bank is produced by the process of claim 3.

10

9. A stable serum-free Vero cell bank having a cell viability of least 80% and a recovery doubling time between 40 and 60 hours after one year, wherein the cell bank is produced by the process of claim 4.

15

ABSTRACT

A serum-free Vero cell banking process provides a standardized and consistent way to generate reliable and stable cell banks for viral vaccine production.

- 5 The substitution of animal-derived substances with plant-derived substances in the growth and freezing media increases the viability of the cells upon thawing and reduces their recovery time, thereby allowing a more precise schedule for manufacturing and more consistent processes.

This Page Is Inserted by IFW Operations
and is not a part of the Official Record

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images may include (but are not limited to):

- BLACK BORDERS
- TEXT CUT OFF AT TOP, BOTTOM OR SIDES
- FADED TEXT
- ILLEGIBLE TEXT
- SKEWED/SLANTED IMAGES
- COLORED PHOTOS
- BLACK OR VERY BLACK AND WHITE DARK PHOTOS
- GRAY SCALE DOCUMENTS

IMAGES ARE BEST AVAILABLE COPY.

**As rescanning documents *will not* correct images,
please do not report the images to the
Image Problem Mailbox.**

FIG. 1: Scale-up scheme for 100 vial cell banks. Cell banks of 200 vials and 400 vials follow the above scheme using two Cell Factories and four Cell Factories respectively.

Sunday	Monday Day 0	Tuesday Day 1	Wednesday Day 2	Thursday Day 3	Friday Day 4	Saturday Day 5
	Thaw two vials of vero cells into a T-150 flask.	Refeed T-150 flask.			Split T-150 flask into five fresh T-150 flasks.	
Day 6	Day 7 Refeed five T-150 flasks.	Day 8	Day 9 Seed Nunc Cell Factory with five T-150 flasks.	Day 10	Day 11 Refeed Nunc Cell Factory.	Day 12
Day 13	Day 14 Refeed Nunc Cell Factory.	Day 15	Day 16 Harvest Nunc Cell Factory. Freeze cell bank.	Day 17	Day 18	Day 19

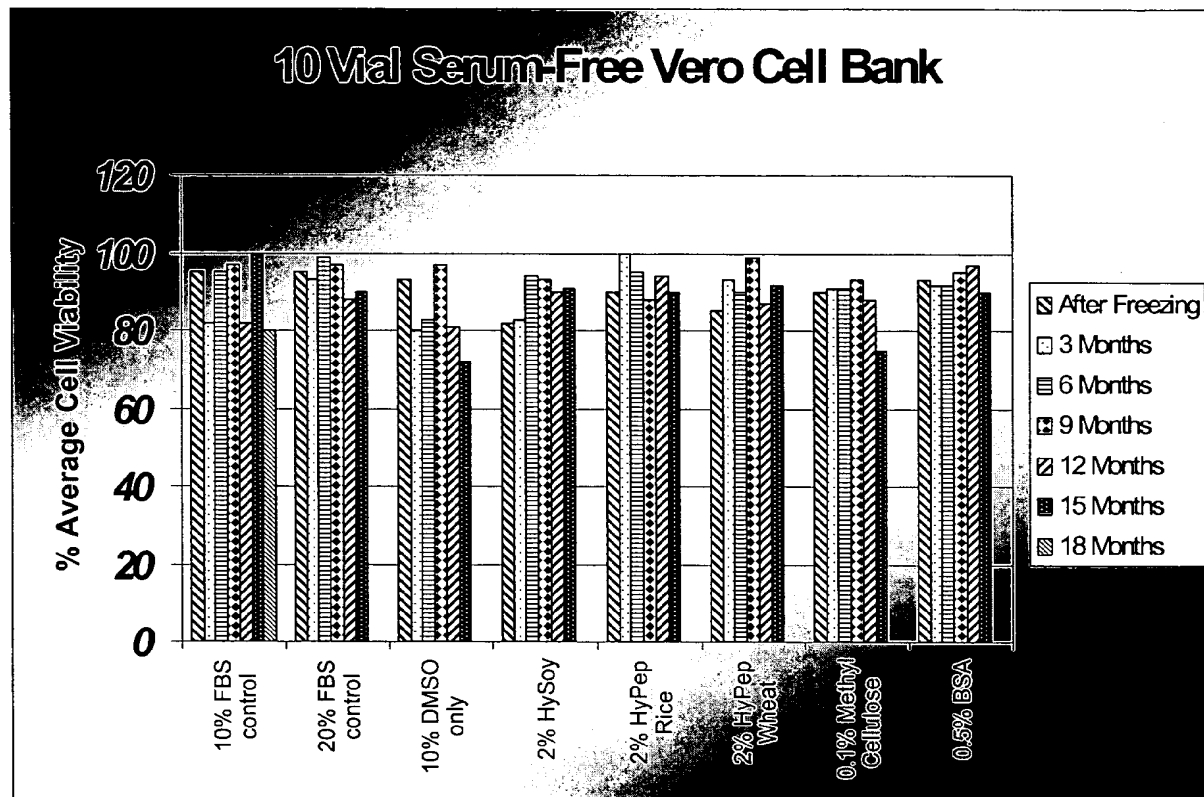


FIG. 3: The average of five samples of each cell bank was taken at various time points. The samples were quickly thawed and placed on T flasks. Cells were counted periodically to determine recovery-doubling times.

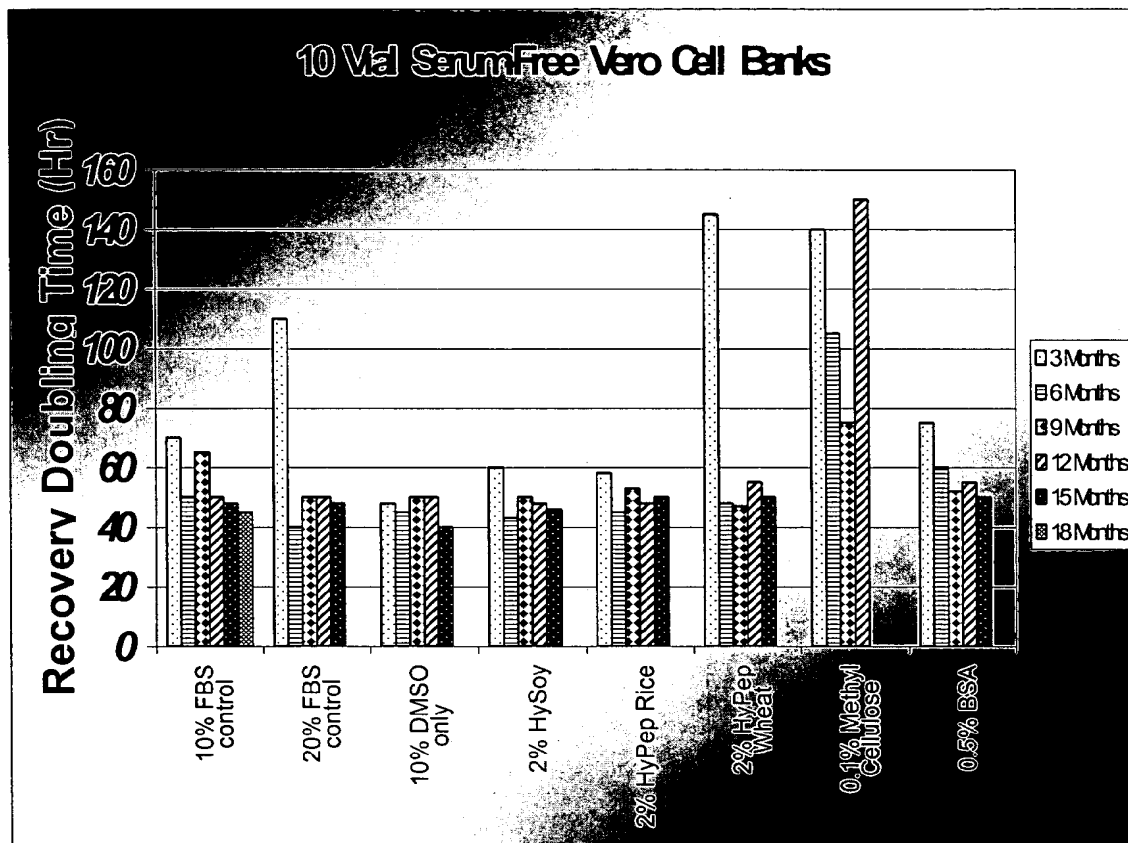


FIG. 4: The average of five samples of each cell bank was taken at various time points. The samples were quickly thawed and placed on T flasks. Cells were counted periodically to determine recovery-doubling times.


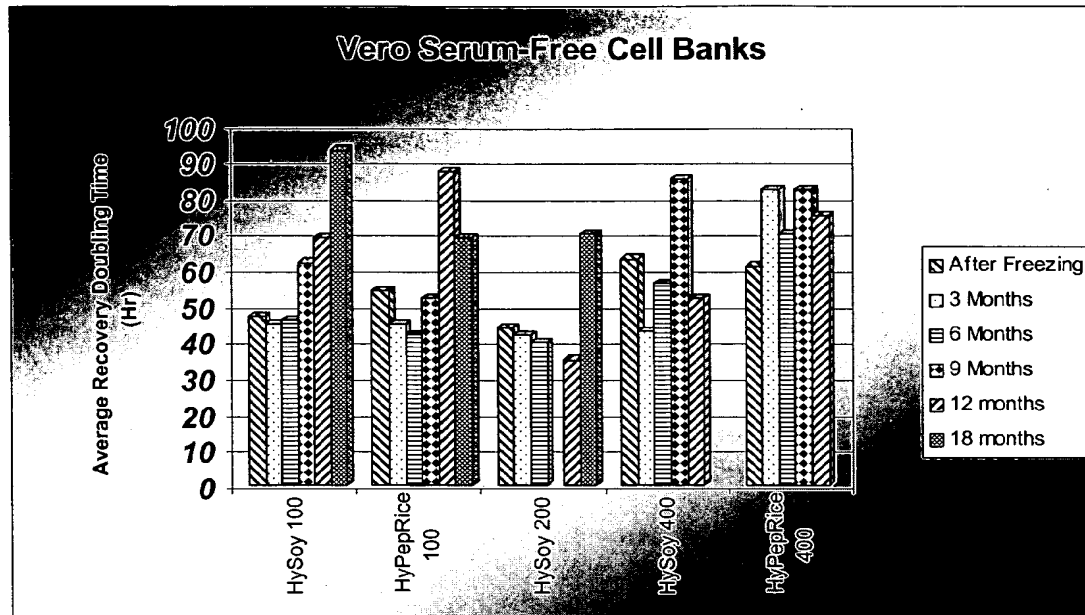


FIG. 5: The average of five samples of each cell bank was taken at various time points. The samples were quickly thawed and viability was determined using Trypan-blue.

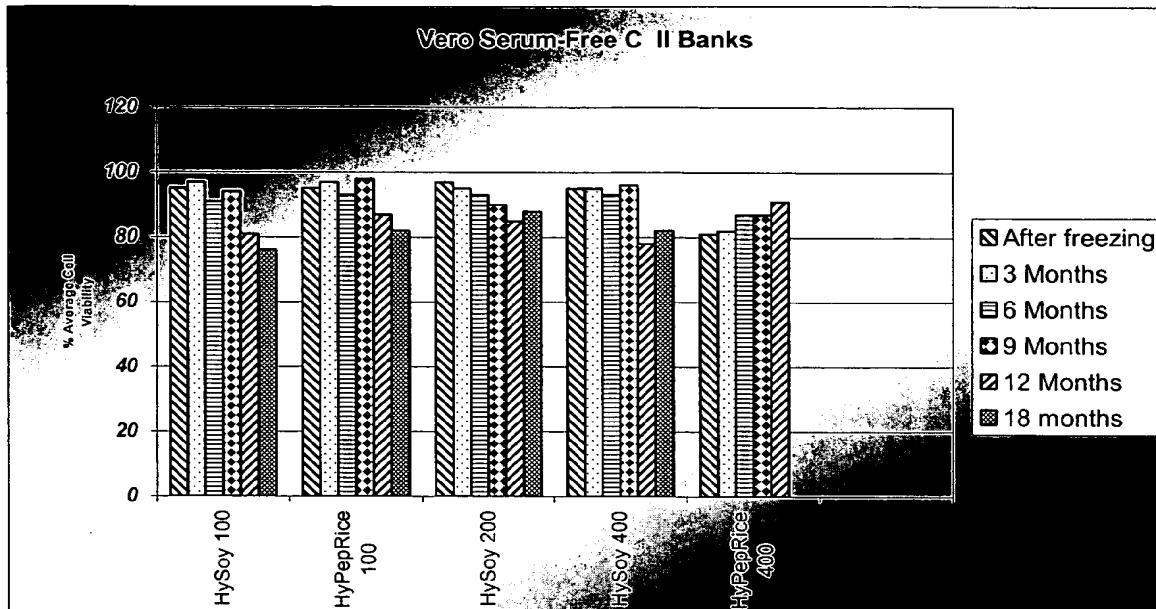

Application Data Sheet

Application Information

Application Type::	Provisional
Subject Matter::	Utility
Suggested Classification::	
Suggested Group Art Unit::	
CD-ROM or CD-R?::	
Number of CD disks::	
Number of copies of CDs::	
Sequence submission?::	
Computer Readable Form (CRF)?::	
Number of copies of CRF::	
Title::	SERUM-FREE VERO CELL BANKING PROCESS
Attorney Docket Number::	AM101242 L1
Request for Early Publication?::	No
Request for Non-Publication?::	No
Suggested Drawing Figure::	Y
Total Drawing Sheets::	5

Applicant Information

Applicant Authority Type::	Inventor
Primary Citizenship Country::	Estonia
Given Name::	Ene
Middle Name::	
Family Name::	Allikmets
City of Residence::	Cornwall on Hudson
State or Province of Residence::	New York
Country of Residence::	USA
Street of mailing address::	200 Maple Road
City of mailing address::	Cornwall on Hudson
State or Province of mailing address::	New York
Country of mailing address::	US
Postal or Zip Code of mailing address::	12520

Correspondence Information

Correspondence Customer Number:: 25291
Phone:: 845-602-4760
Fax Number:: 973-660-7973
E-Mail address:: RozekC@Wyeth.com

Representative Information

Representative Customer Number::	25291
----------------------------------	-------

Domestic Priority Information

Application::	Continuity Type::	Parent Application::	Parent Filing Date::
This Application			

Foreign Priority Information

Country::	Application Number::	Filing Date::	Priority Claimed::

Assignee Information

Assignee Name:: Wyeth
Street of mailing address:: Five Giralda Farms
City of mailing address:: Madison
State or Province of mailing address:: New Jersey
Country of mailing address:: USA
Postal or Zip Code of mailing address:: 07940